


Kontakt™ W

Simple and reliable
to increase efficiency

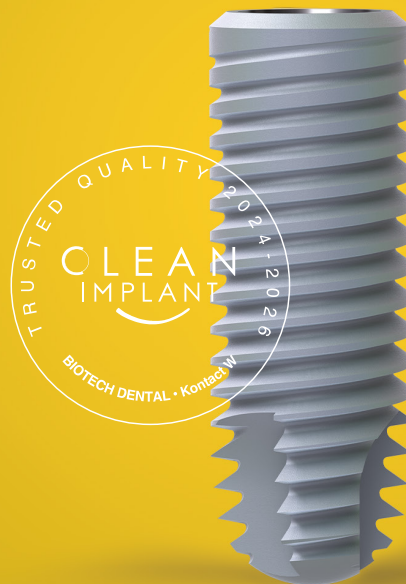

BIOTECH DENTAL

MORE THAN 35 YEARS OF EXPERTISE

Designed with scientific knowledge in implantology established for more than 35 years, and current manufacturing requirements, the Kontakt™ W implant increase technical performances with a high quality level.

Its purpose is to facilitate implantology access, both for practitioners and patients.



✓ **Narrowed, chamfered and micro-structured implant neck**

The narrowing neck reduces bone compression of the cortical bone. The micro-structured chamfer promotes the retention of blood coagulum which is necessary for bone reconstruction.




✓ **Platform switching**

Platform switching is an important factor concerning the tissue stability. The aim is to limit the occurrence of peri-implantitis.



✓ **Conical connection with hexagonal indexing**

Improved bacterial seal and eliminates micromovements between the implant and the abutment. The implant-abutment mechanical strength is enhanced. Placement of prosthetic parts becomes easier.

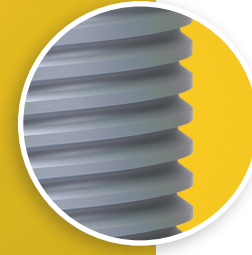

✓ **Performing implant surface with grade IV titanium (T60)**

The etched sandblasted surface improve wettability and osseointegration. The Kontakt™ W implant benefits from the same surface as the Kontakt™ S implant, which is certified CleanImplant.



✓ **Atraumatic spherical apex**

Thanks to its design, the apex protects the anatomical structure during surgery and allows a less invasive osteotomy.

2

TREATMENT SIMPLICITY AND FLEXIBILITY

A SINGLE CONNECTION
FOR THE WHOLE RANGE
OF KONTACT™ W IMPLANT


3 diameters and 5 different
lengths.

> **3 body diameters:**



> **5 lengths:**

6 mm / 8 mm / 10 mm / 12 mm / 14 mm

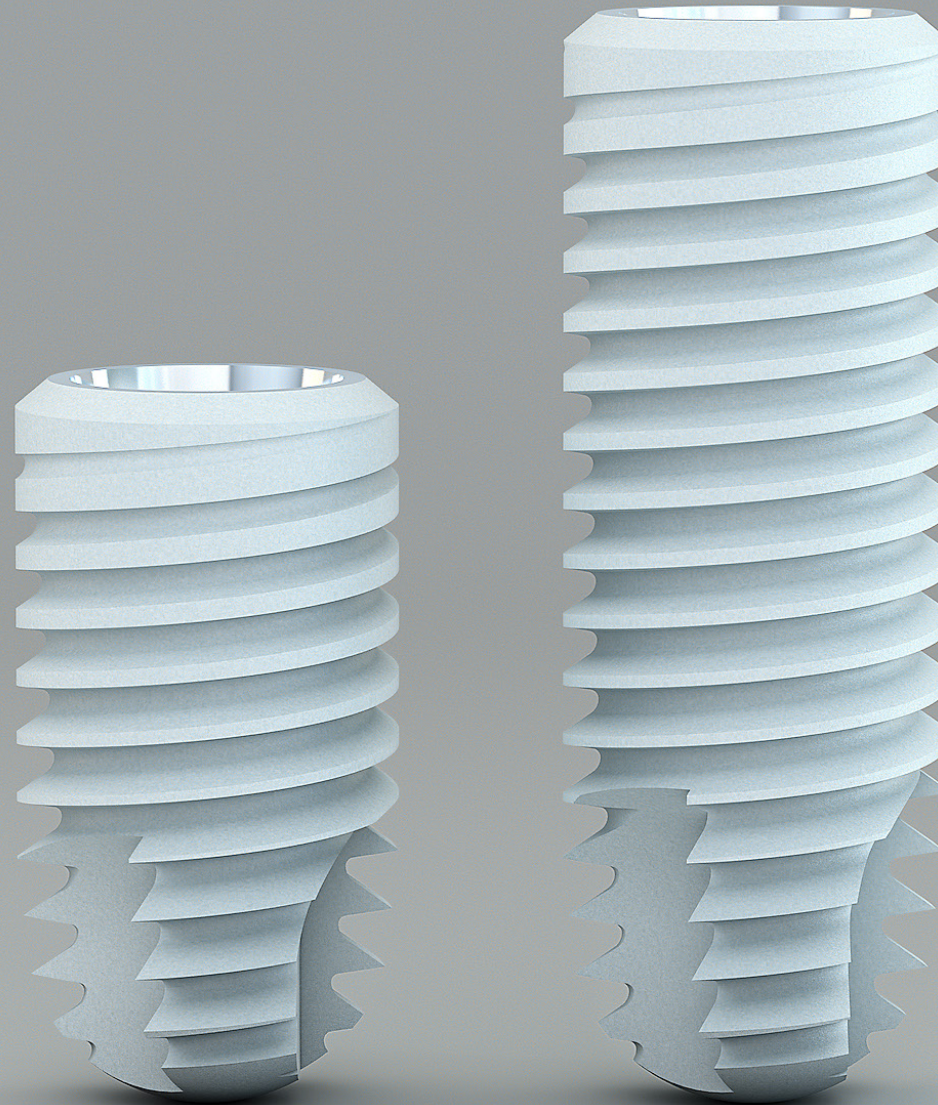


INCLUDE IMPLANT
SCANNING INTO YOUR
DIGITAL WORKFLOW


- ✓ Improve your digital workflow, your patients communication and accelerate the planning of your cases.
- ✓ Make your impressions easier with an intraoral scanner.

3

A SINGLE SURGICAL KIT FOR ALL KONTACT™ IMPLANTS

- ✓ **Compact** for optimal storage.
- ✓ **Understandable** thanks to colour markings for quick identification of instruments.
- ✓ **Practical** with fast and easy opening for instruments access.
- ✓ **Complete disassembly** for thorough cleaning and sterilisation.

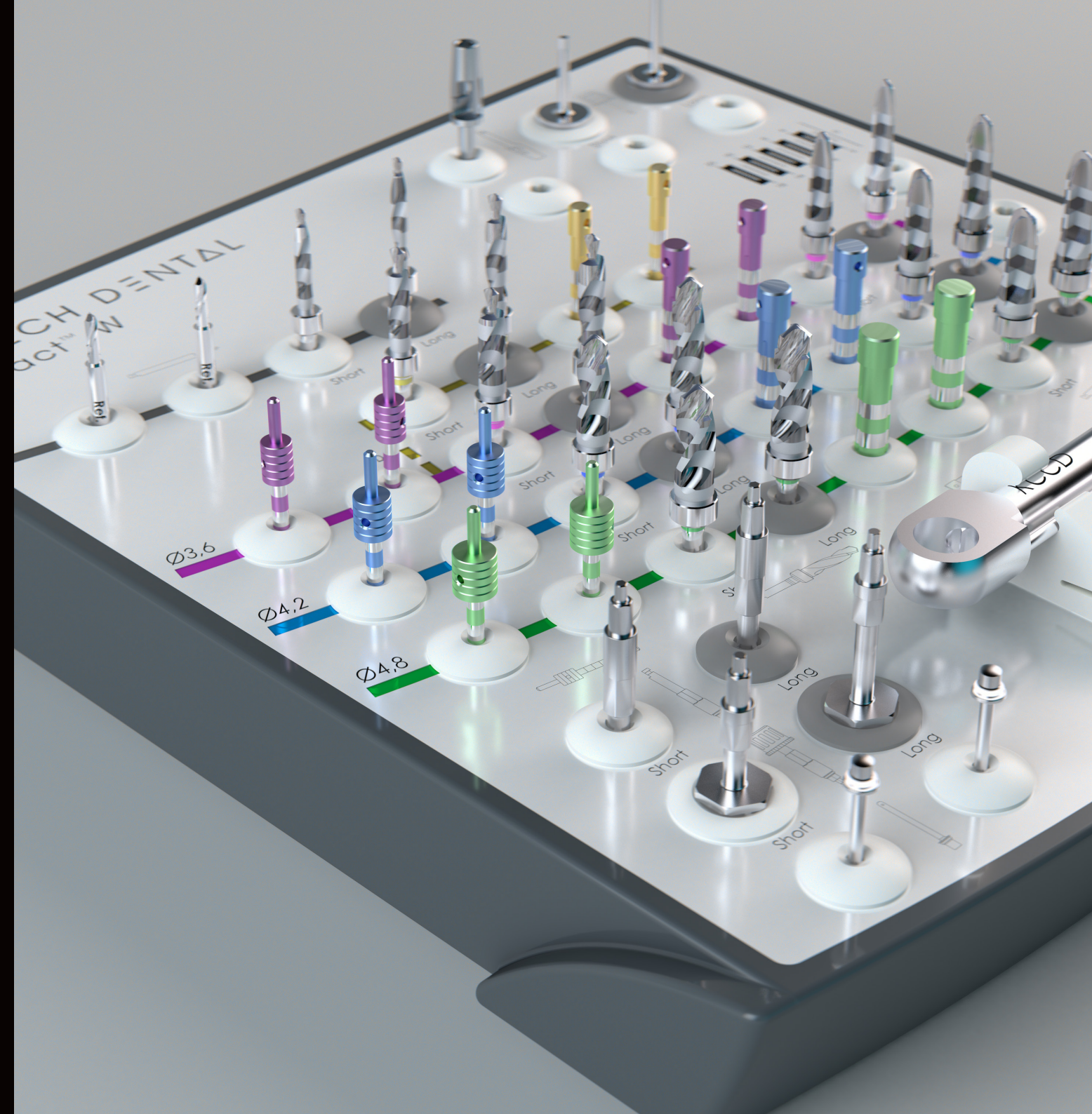

Lengths	Diameters		
	Ø 3.6 mm	Ø 4.2 mm	Ø 4.8 mm
6mm		 W4206	 W4806
8mm	 W3608	 W4208	 W4808
10mm	 W3610	 W4210	 W4810
12mm	 W3612	 W4212	 W4812
14mm	 W3614	 W4214	 W4814

4

A WIDE RANGE OF PROSTHETIC SOLUTIONS

- ✓ Cemented prosthesis
- ✓ Screw-retained prosthesis
- ✓ CAD-CAM prosthesis



NOTES

NOTES



BIOTECH DENTAL



305, Allées de Craponne
13300 Salon-de-Provence - France
Tel.: +33 (0)4 90 44 60 60
Fax: +33 (0)4 90 44 60 61

CustomerCareService@biotech-dental.com
www.biotech-dental.com

www.biotechgalaxy.com



AFNOR Cert. 73017

Biotech Dental - 305, Allées de Craponne - 13300 Salon-de-Provence, France.

Biotech Dental - A simplified company limited by shares (S.A.S.) with a capital of 24,866,417€. Corporate and Trade Register of Salon de Provence: 795 001 304 - SIRET: 795 001 304 00018 - VAT N°: FR 31 79 500 13 04.

Kontakt™ W: manufactured by Biotech Dental. Class I, IIa and IIb medical devices dedicated to dental implantology. CE0477.

All of these medical devices must be used by qualified and trained health professionals. All brands are under the responsibility of their respective manufacturers. Read the instructions carefully. No refund from French social security. Non contractual visuals. Do not litter in public areas.

Printed by VALLIERE - 163, Avenue du Luxembourg - ZAC des Molières - 12140 Miramas - France.